


# Meet the next member of the family.

## The Baum Desk Aneroid

### FEATURES:

- 300 mm Hg Aneroid instrument mounted at a 20° angle for ease of viewing.
- Steel stand finished with Neutral Gray textured baked enamel and non-skid base.
- Complies with ANSI/AAMI SP-10 for accuracy and performance (1% full scale  $\pm$  3 mmHg).
- Neutral White ABS plastic manometer housing with chrome plated steel bracket.
- Unique high contrast scale and red pointer to increase visibility and reduce parallax error.
- Manufactured with a low-friction precision movement and twin pre-seasoned phosphor-bronze diaphragm capsules.

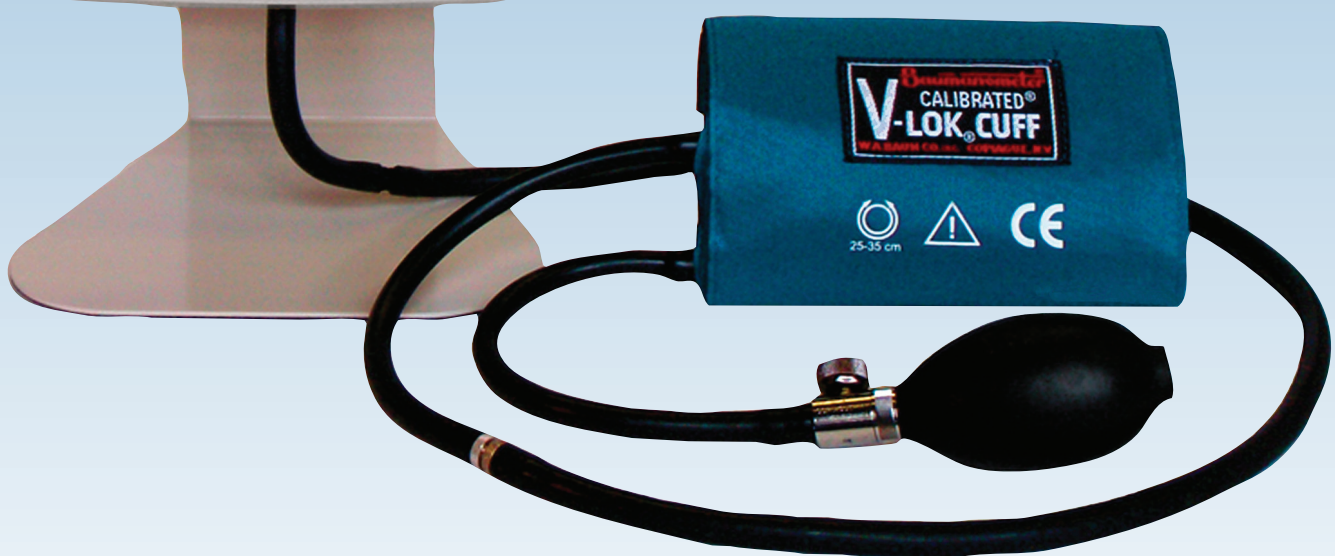

**W.A. Baum Co. Inc.**

Since 1916

Originator and Maker of Accurate Blood Pressure Apparatus Exclusively

620 Oak Street, Copiague, NY 11726 U.S.A.

www.wabaum.com • info@wabaum.com • Toll-Free: 1-888-281-6061 • Tel: 631-226-3940 • Fax: 631-226-3969

► Indicates a frequently ordered product

# Features and Specifications



1820

## Calibrated® V-Lok® Inflation System:

**Consists of:** Bulb, Air-Flo® Control Valve, Inflation Bag and Calibrated® V-Lok® Cuff.

**Accurate:** The Calibrated® V-Lok® System immediately indicates whether the correct size cuff is being used. The center of the bladder is clearly marked on the inner and outer surface of the cuff for correct application on either the right or left limb.


**Comfortable:** Cuffs conform to the shape of the limb and have no hard, stiff edges.

**Durable:** Baum cuffs are made of tightly woven, heat set Dacron™ polyester fabric. They are double stitched at critical points for long wear and treated with an effective anti-microbial agent. Hook and loop fasteners of matched holding strength provide the optimal number of open-close cycles.

**Non-Latex:** Available for those who are hypersensitive to natural rubber latex.

**Guaranteed:** Latex and **Non-Latex** inflation bags, bulbs and tubing are guaranteed for five years.

For additional information, refer to the Calibrated® V-Lok® Cuff Data Sheet or visit [www.wabaum.com](http://www.wabaum.com).

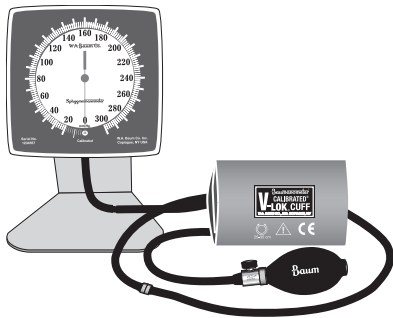

0920

## Baum Desk Aneroid:

- 300 mm Hg Aneroid instrument mounted at a 20° angle for ease of viewing.
- Complies with ANSI/AAMI SP-10 for accuracy and performance (1% full scale ± 3 mm Hg).
- Neutral White ABS plastic manometer with chrome plated steel bracket.
- Unique high contrast scale and red pointer increase visibility and reduce parallax error.
- Manufactured with a low-friction precision movement and twin pre-seasoned phosphor-bronze diaphragm capsules.

**Steel Stand Finish:** Neutral Gray textured baked enamel with non-skid base.

**Packed Dimensions:** 10-5/8" L x 7-1/4" W x 7-1/2" D. **Packed Weight:** 3-1/4 lbs.

## Desk Aneroid Instrument Ordering Information:

**Includes:** Calibrated® V-Lok® System with 25" of straight tubing and a steel stand.

**Cat. No.**

| Latex  | Non-Latex | Inflation System                     | Instrument Size          |
|--------|-----------|--------------------------------------|--------------------------|
| ► 0920 | 0920NL    | Adult Calibrated® V-Lok®             | 6-1/4" H x 6-1/4"W x 4"D |
| 0921   | 0921NL    | Child/Small Adult Calibrated® V-Lok® | 6-1/4" H x 6-1/4"W x 4"D |
| 0922   | 0922NL    | Infant Calibrated® V-Lok®            | 6-1/4" H x 6-1/4"W x 4"D |
| 0923   | 0923NL    | Newborn Calibrated® V-Lok®           | 6-1/4" H x 6-1/4"W x 4"D |
| 0924   | 0924NL    | Thigh Calibrated® V-Lok®             | 6-1/4" H x 6-1/4"W x 4"D |
| 0925   | 0925NL    | Large Arm Calibrated® V-Lok®         | 6-1/4" H x 6-1/4"W x 4"D |

## Guarantee:

Every Baum Aneroid instrument is guaranteed against defects in workmanship and material for a period of one year from date of purchase. This guarantee is invalid when other than genuine W.A. Baum parts are used. Latex and latex-free components are guaranteed for 5 years.

We recommend that every aneroid instrument be checked for accuracy periodically by following our "Testing Instrument Accuracy" procedure\* to determine if it needs adjustment. Adjustments will be made to Baum aneroid instruments without charge when returned, postage prepaid, unless the instrument has been damaged in shipment, tampered with or abused. Our pledge is that our customers will be completely satisfied with our products.

\*Printed copies of "Testing Instrument Accuracy" are available on request or by visiting the FAQ section of our website at: [www.wabaum.com](http://www.wabaum.com).

Baumanometer, Standby, V-Lok, Extendex, Calibrated, Hi/Lo, Air-Flo and Roll-By are registered trademarks of W.A. Baum Co., Inc.

**W.A. Baum Co., Inc.**  
620 Oak Street  
Copiague, NY 11726 USA

Tel: 631-226-3940  
Fax: 631-226-3969  
Toll Free: 1-888-281-6061

Web: [www.wabaum.com](http://www.wabaum.com)  
Email: [info@wabaum.com](mailto:info@wabaum.com)

**Since 1916 Originator and Maker of Accurate Blood Pressure Apparatus Exclusively**

# W.A. Baum Co. Inc.