

The Only Pocket Aneroid That Measures Up To Our Standards.


W.A. Baum Company has specialized in the manufacture of scientifically accurate clinical sphygmomanometers for more than 85 years. Each instrument is built to agree with our master



Baumanometer® Blood Pressure Standard The World Over.

instruments, which are traceable to the United States National Institute for Standards and Technology.

The Baum Pocket Aneroid complies with ANSI/AAMI SP-10 for accuracy and performance (1% full scale \pm 3 mmHg). It is equipped with the World Standard Calibrated® V-Lok® Cuff and our exclusive Bulb and Air-Flo® Control Valve. For more information on the Baum Pocket Aneroid call Toll-Free 1-888-281-6061 or e-mail at info@wabaum.com for assistance in locating an authorized Baum dealer.

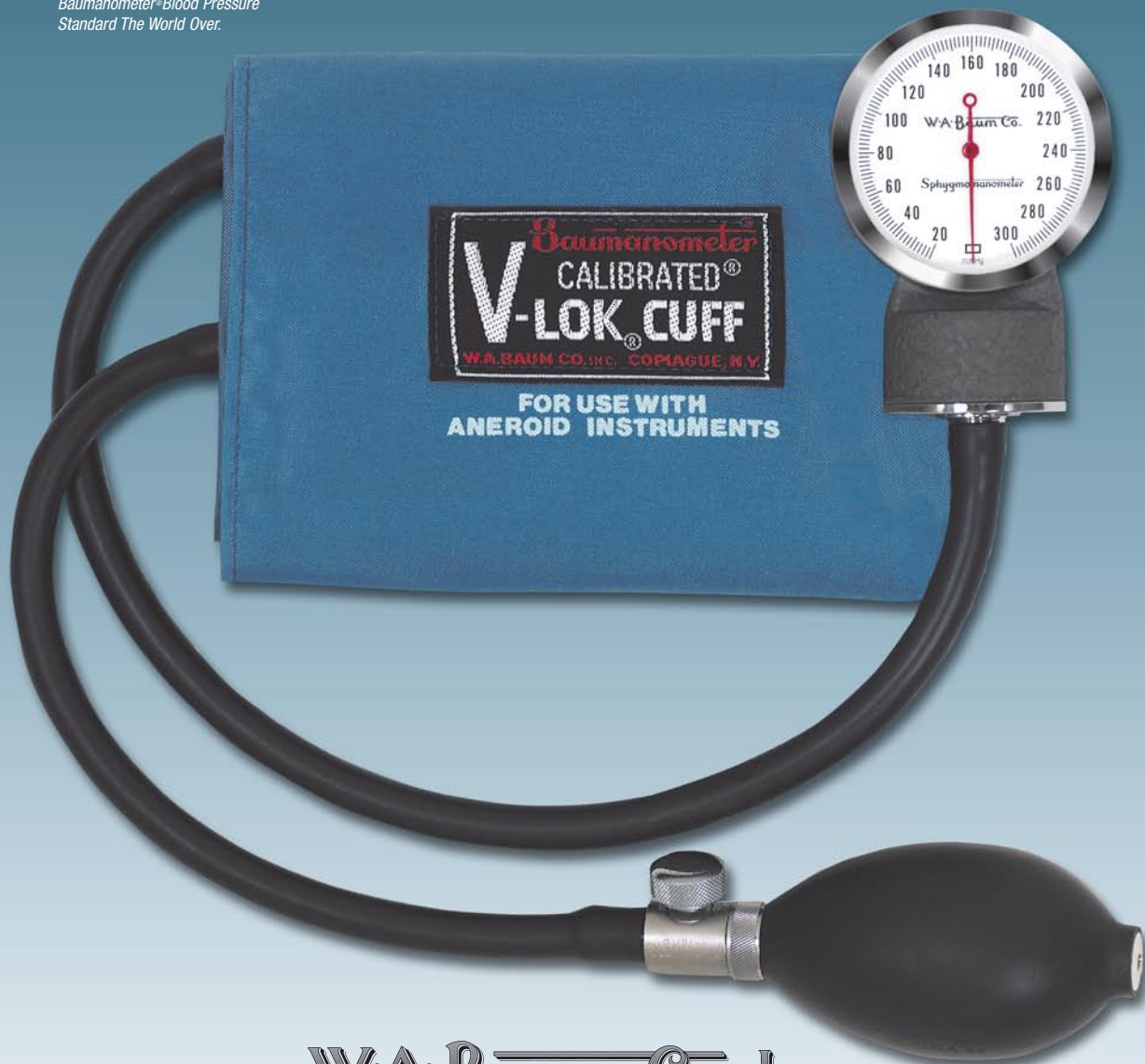

W.A. Baum Co. Inc.

SINCE 1916


620 Oak Street, Copiague, NY 11726 USA

www.wabaum.com • info@wabaum.com • Toll-Free: 1-888-281-6061 • Tel: 631-226-3940 • Fax: 631-226-3969

Features and Specifications



1820AC

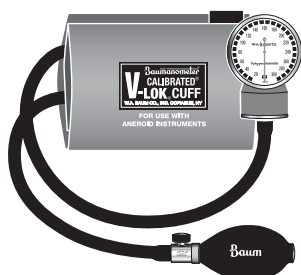

1050



Zippered Case

Calibrated® V-Lok® Inflation System:

Accurate: The Calibrated® V-Lok® system clearly indicates whether or not the correct size cuff is being used. The center of the bladder is clearly marked on the inner and outer surfaces of the cuff for correct application on either the right or left limb.

Comfortable: Cuffs conform to the shape of the limb and have no hard, stiff edges.

Durable: Baum cuffs are made of tightly woven, heat set Dacron™ polyester fabric. They are double stitched at critical points for long wear and treated with an effective anti-microbial agent.

Non-Latex: Available for those who are hypersensitive to natural rubber latex.

Guaranteed: Latex and **Non-Latex** inflation bags, bulbs and tubing are guaranteed for five years.

Pocket Aneroid:

The Baum Pocket Aneroid is a 300 mmHg portable instrument. It's built to agree with our master instruments, which are traceable to the United States National Institute of Standards and Technology. It complies with ANSI/AAMI SP-10 for accuracy and performance (1% full scale ± 3 mmHg). The unique, high contrast scale and red pointer increase visibility and reduce parallax error. They are manufactured with a precision low friction movement and twin pre-seasoned phosphor-bronze diaphragm capsules.

Materials and Finish:

The die cast zinc manometer housing is finished in warm gray textured baked enamel. The standard Adult Calibrated® V-Lok® Cuff is medium blue Dacron® polyester fabric. The adult size is also available in dark brown, dark green, scarlet, tan and gray.

Pocket Aneroid Ordering Information:

Cat. No.

Latex	Non-Latex	Model	Inflation System	Calibration	Holder
1050	1050NL	Pocket	Adult Calibrated® V-Lok®	300 mm Hg	Zippered Case
1051	1051NL	Pocket	Child Calibrated® V-Lok®	300 mm Hg	Zippered Case
1052	1052NL	Pocket	Infant Calibrated® V-Lok®	300 mm Hg	Zippered Case
1054	1054NL	Pocket	Thigh Calibrated® V-Lok®	300 mm Hg	None
1055	1055NL	Pocket	Large Arm Calibrated® V-Lok®	300 mm Hg	Zippered Case

Guarantee:

Every Baum Aneroid instrument is guaranteed against defects in workmanship and material for a period of one year from date of purchase. This guarantee is invalid when other than genuine W.A. Baum parts are used. Latex and latex-free components are guaranteed for 5 years.

We recommend that every aneroid instrument be checked for accuracy periodically by following our “[Testing Instrument Accuracy](#)” procedure* to determine if it needs adjustment. Adjustments will be made to Baum aneroid instruments without charge when returned, postage prepaid, unless the instrument has been damaged in shipment, tampered with or abused. Our pledge is that our customers will be completely satisfied with our products.

*Printed copies of “[Testing Instrument Accuracy](#)” are available on request or by visiting the FAQ section of our website at: www.wabaum.com.

Baumanometer, Standby, V-Lok, Extendex, Calibrated, Hi/Lo, Air-Flo and Roll-By are registered trademarks of W.A. Baum Co., Inc.

W.A. Baum Co., Inc.

620 Oak Street
Copliague, NY 11726 USA

Tel: 631-226-3940

Fax: 631-226-3969

Toll Free: 1-888-281-6061

Web: www.wabaum.com

Email: info@wabaum.com

Since 1916 Originator and Maker of Accurate Blood Pressure Apparatus Exclusively

W.A. Baum Co. Inc.