



Wall Aneroid

Looking for an alternative?

Try the **Wall Aneroid** from W.A. Baum. It's built to agree with our master instruments, which are traceable to the United States National Institute of Standards and Technology. It complies with AAMI SP-10 for accuracy and performance (1% full scale \pm 3 mmHg). Its unique, high contrast scale and red pointer increase visibility and reduce parallax error.

The manometer is manufactured with a metal mounting bracket, high impact ABS plastic manometer housing, precision low-friction movement and twin pre-seasoned phosphor-bronze diaphragm capsules.

Wall Aneroids can be ordered with any of our Calibrated® V-Lok® Cuffs, 8 feet of Extendex® coiled tubing, nylon coated steel inflation system storage basket and superior fasteners for wall mounting.

For additional information on the Baum **Wall Aneroid** or for assistance in locating an authorized Baum dealer contact us at: info@wabaum.com

W.A. Baum Co. Inc.

SINCE 1916


620 Oak Street, Copiague, NY 11726 USA
www.wabaum.com • info@wabaum.com

Toll-Free: 1-888-281-6061 • Tel: 631-226-3940 • Fax: 631-226-3969


► Indicates a frequently ordered product.



1820

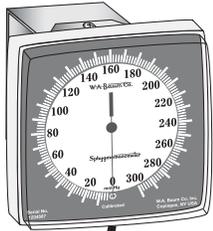

0950

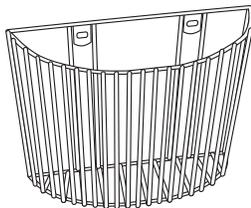

2422

Specifications

Calibrated® V-Lok® Inflation System:

Accurate: The Calibrated® V-Lok® system clearly indicates whether or not the correct size cuff is being used. The center of the bladder is clearly marked on the inner and outer surfaces of the cuff for correct application on either the right or left limb.

Comfortable: Cuffs conform to the shape of the limb and have no hard, stiff edges.

Durable: Baum cuffs are made of tightly woven, heat set Dacron™ polyester fabric. They are double stitched at critical points for long wear and treated with an effective anti-microbial agent.

Non-Latex: Available for those who are hypersensitive to natural rubber latex.

Guaranteed: Latex and **Non-Latex** inflation bags, bulbs and tubing are guaranteed for five years.

Wall Aneroid:

The Baum Wall Aneroid is a 300 mmHg swivel-type instrument (110° arc) for wall mounting. The manometer housing is neutral white ABS plastic. The wall mounting bracket is chrome plated steel and together they measure 6.25"H x 6.25"W x 4"D.

The standard Adult Calibrated® V-Lok® Cuff is medium blue Dacron® polyester. The adult size is also available in dark brown, dark green, scarlet, tan and gray.

The inflation system storage basket is neutral white nylon coated steel and measures 7"H x 6.5"W x 3.75"D. The large storage basket is 7"H x 9.75"W x 3.75"D.

Complete instruments include: Calibrated® V-Lok® Inflation System, nylon coated steel storage basket, 8 feet of Extendex® coiled tubing and superior fasteners for wall mounting.

Ordering Information:

Cat. No. Description

- **0950** Wall Aneroid w/Adult Calibrated® V-Lok® and basket
- 0950NL** Wall Aneroid w/Adult Non-Latex Calibrated® V-Lok® and basket
- 0951** Wall Aneroid w/Child/Small Adult Calibrated® V-Lok® and basket
- 0952** Wall Aneroid w/Infant Calibrated® V-Lok® and basket
- 0955** Wall Aneroid w/Large Arm Calibrated® V-Lok® and basket
- 0960** Wall Aneroid w/Adult Calibrated® V-Lok® (**no basket**)

Note:

To order a Wall Aneroid with a size or type of inflation system that is not listed, use **Cat. No. 0933**, select the separate catalog number for the size or type of inflation system needed from **page 8** in our catalog and your choice of storage basket. Each item will be listed separately on your order. Order Non-Latex items by adding NL to the catalog number.

Cat. No. Description

- 0933** Wall Aneroid with 8 ft. Extendex® tubing only
- 0970** Wall Aneroid w/8 ft. Extendex® tubing, large storage basket (#2422) and three complete Calibrated® V-Lok® Inflation Systems (Adult, Child/Small Adult and Large Arm) packed in one carton.
- 2420** Storage basket (nylon coated steel)
- 2422** Large storage basket (nylon coated steel)

Guarantee:

Every Baum Aneroid instrument is guaranteed against defects in workmanship and material for a period of one year from date of purchase. This guarantee is invalid when other than genuine W.A. Baum parts are used.

We recommend that every aneroid instrument be checked for accuracy periodically by following our "Testing Instrument Accuracy" procedure* to determine if it needs adjustment. Adjustments will be made to Baum aneroid instruments without charge when returned, postage prepaid, unless the instrument has been damaged in shipment, tampered with or abused. Our pledge is that our customers will be completely satisfied with our products.

*Printed copies of "Testing Instrument Accuracy" are available on request or by visiting the FAQ section of our website at: www.wabaum.com.

Baumanometer, Standby, V-Lok, Extendex, Calibrated, Hi/Lo, Air-Flo and Roll-By are registered trademarks of W.A. Baum Co., Inc.

W.A. Baum Co., Inc.

620 Oak Street
Copliague, NY 11726 USA

Tel: 631-226-3940

Fax: 631-226-3969

Toll Free: 1-888-281-6061

Web: www.wabaum.com

Email: info@wabaum.com

Since 1916 Originator and Maker of Blood Pressure Apparatus Exclusively

W.A. Baum Co. Inc.