
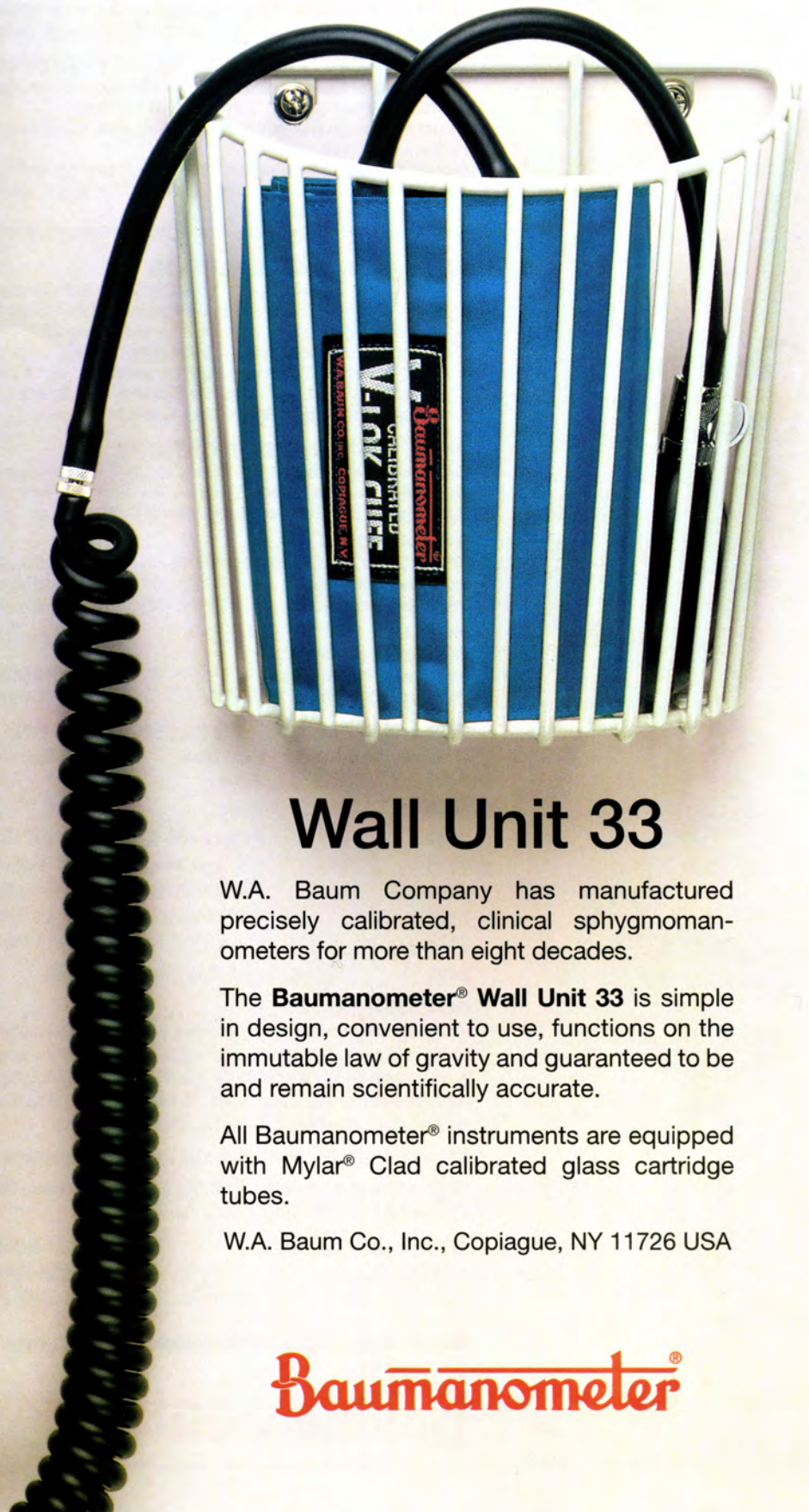


# World Standard

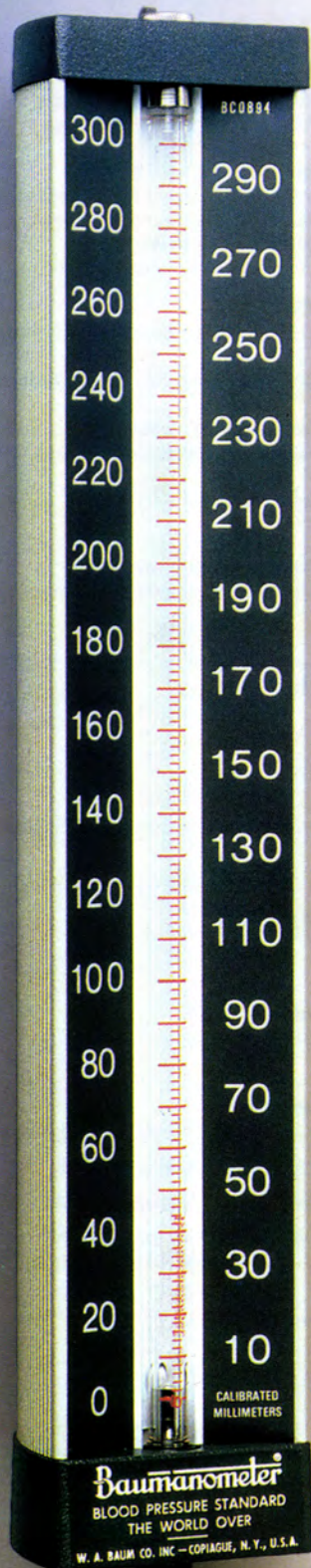

## Wall Unit 33

W.A. Baum Company has manufactured precisely calibrated, clinical sphygmomanometers for more than eight decades.

The **Baumanometer® Wall Unit 33** is simple in design, convenient to use, functions on the immutable law of gravity and guaranteed to be and remain scientifically accurate.


All Baumanometer® instruments are equipped with Mylar® Clad calibrated glass cartridge tubes.

W.A. Baum Co., Inc., Copiague, NY 11726 USA


**Baumanometer®**




# Features & Specifications

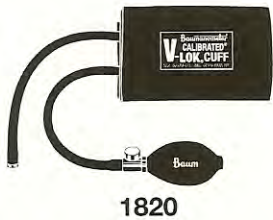

1820

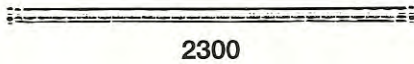


2300

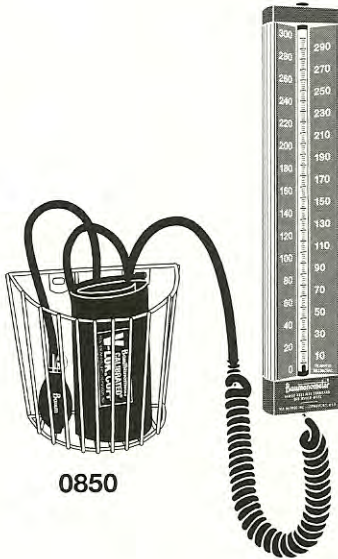

0850



2420

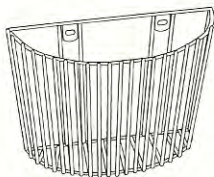

2422



2421

## CALIBRATED® V-LOK® INFLATION SYSTEM

**Accurate:** The Calibrated® V-Lok® system clearly indicates whether or not the correct size cuff is being used. The center of the bladder is clearly marked on the inner and outer surfaces of the cuff for correct application on either the right or left limb.

**Comfortable:** Cuffs conform to the shape of the limb and have no hard, stiff edges.

**Durable:** Baum cuffs are made of tightly woven, heat set Dacron™ polyester fabric. They are double stitched at critical points for long wear and treated with an effective anti-microbial agent.

**Non-Latex:** Available for those who are hypersensitive to natural rubber latex.

**Guaranteed:** Latex and Non-Latex inflation bags, bulbs and tubing are guaranteed for five years.

## MYLAR® CLAD CALIBRATED GLASS CARTRIDGE TUBE

**Safe:** Several layers of clear Mylar® film have been applied to strengthen the tube and maintain its structural integrity even if the glass is broken.

**Effective:** Large bore tube (5 mm) prevents separation of mercury and is 25% larger than required by Federal Specifications.

**Durable:** Cartridge tube end faces are flame-glazed to prevent chipping; ceramic graduation marks are fused to the glass for permanence. The entire tube is carefully annealed to eliminate both longitudinal and ring strain.

**Guaranteed:** Free replacement of the cartridge tube if it is ever broken.

## WALL UNIT 33 DESCRIPTION

- 300 mmHg swivel-type (110° arc) clinical sphygmomanometer for wall mounting.
- Meets the American National Standard ANSI/AAMI SP-10 for accuracy and performance.
- Available with Latex or Non-Latex components.
- Extruded aluminum alloy manometer housing with diecast zinc top and base.
- 3 convenient cuff storage options available: large basket, small basket or stainless steel wall clip.
- Calibrated® V-Lok® Inflation System with 8 feet of Extendex® Tubing.
- Mylar® Clad Calibrated Glass Cartridge Tube.

## WALL UNIT 33 SPECIFICATIONS

**Instrument Finish:** Neutral White textured baked enamel (standard)

**Optional Instrument Colors:** Silvertone, Bright White, Pastel Blue, Pastel Rose and Pastel Green

**Cuff Color:** Medium Blue (standard) **Optional Cuff Colors:** Dark Brown, Dark Green, Gray, Scarlet and Tan (Adult Size Only)

### WITH BASKET:

**Basket Color:** Neutral White (standard)

**Instrument Dimensions:** 14 3/4" L x 2 1/2" W x 3" H

**Packaged Dimensions:** 17" L x 7" W x 4" H

**Optional Basket Colors:** Chrome and Bright White

**Instrument Weight:** 2.25 lbs.

**Packaged Weight:** (w/small basket) 4.20 lbs.

### WITH CLIP:

**Wall Clip Finish:** Stainless Steel

**Instrument Dimensions:** 14 3/4" L x 2 1/2" W x 3" H

**Packaged Dimensions:** 16" L x 4 1/4" W x 4 3/4" H

**Instrument Weight:** 2.25 lbs.

**Packaged Weight:** 3.45 lbs.

## WALL UNIT 33 ORDERING INFORMATION:

Cat. No.	Non-Latex	Cuff Size	Inflation System	Holder
● 0850	0850NL	Adult	Calibrated® V-Lok®	Wall Basket (Nylon Coated Steel)
0820	0820NL	Adult	Calibrated® V-Lok®	Wall Clip (Stainless Steel)
0833	0833NL	<b>None (with 8 ft. Extendex® Tubing Only). No Inflation System. No Holder.</b>		
● 0870	0870NL	<b>Wall Unit 33 w/Large Wall Basket (#2422) and three complete Calibrated® V-Lok® Inflation Systems (Adult, Child/Small Adult and Large Arm) packed in one carton.</b>		
2420		Wall Basket (Nylon Coated Steel)		6 1/2" W x 7" H x 3 3/4" D
2422		Large Wall Basket (Nylon Coated Steel)		9 3/4" W x 7" H x 3 3/4" D
2421		Wall Clip (Stainless Steel)		

## GUARANTEES

Every Baumanometer® instrument is a precisely calibrated true gravity manometer. It is guaranteed to be and remain scientifically accurate. The entire apparatus is guaranteed against defects in workmanship and material for a period of one year from date of purchase. We guarantee the free replacement of the Calibrated Cartridge Tube if it is ever broken. Latex and Non-Latex inflation bags, bulbs and tubing are guaranteed for five years.

These guarantees are invalid when other than genuine Baumanometer® parts are used.

Baumanometer, Standby, Roll-By, V-Lok, Extendex, Calibrated, Hi/Lo, Air-Flo are registered trademarks of W.A. Baum Co., Inc.

**W.A. Baum Co., Inc.**

620 Oak Street  
Copiague, NY 11726 USA

Tel: 631-226-3940

Fax: 631-226-3969

Toll Free: 1-888-281-6061

Web: [www.wabaum.com](http://www.wabaum.com)

Email: [info@wabaum.com](mailto:info@wabaum.com)

**Since 1916 Originator and Maker of Accurate Blood Pressure Apparatus Exclusively**

# Baumanometer®