

Guarantee:

Baum Dacron® polyester cuffs are guaranteed against defects in workmanship and material for a period of one year.

Baum Inflation Bags, Bulbs & Tubing are guaranteed for five years.



W.A. Baum Co., Inc.
620 Oak Street
Copiague, NY 11726
U.S.A.



Hawksley & Sons
Marlborough Road
Lancing, West Sussex BN15 8TN
UK
+44 (0) 1903 7528 15
enquiries@Hawksley.co.uk

CE 0197

± 3 mmHg

EN 1060-1:1995
EN 1060-2:1995

Available from your distributor

Bariatric patients have special needs.



Baum Bariatric Inflation System

W.A. Baum Co., Inc.
620 Oak Street
Copiague, NY 11726
www.wabaum.com

Toll Free: 1- 888-281-6061
Phone: 631-226-3940
Facsimile: 631-226-3969
Email: info@wabaum.com

W.A. Baum Co. Inc.

SINCE 1916

ORIGINATOR AND MAKER OF ACCURATE
BLOOD PRESSURE APPARATUS EXCLUSIVELY

The Baum Bariatric Inflation System

Bariatric patients have special needs.

Often the upper arm of a bariatric patient is conical or V-shaped. This inhibits the smooth application of a standard (rectangular) blood pressure cuff because they do not conform to the conical shape of the limb. A gap at the bottom of the cuff near the elbow is a common occurrence and is a clear indication of a misapplied cuff.

THE CHALLENGE:

Correct cuff application on either right or left limb and even compression of the artery for Bariatric patients.

THE SOLUTION:

The Baum Bariatric Inflation System:



- Easy to Apply
- Conforms Naturally to a Conical Limb.

Cat. No. Description

1830	Bariatric Inflation System
1830NL	Bariatric Non-Latex Inflation System

Range of Limb Circumferences 

46 - 66 cm

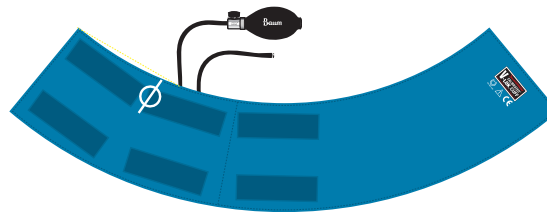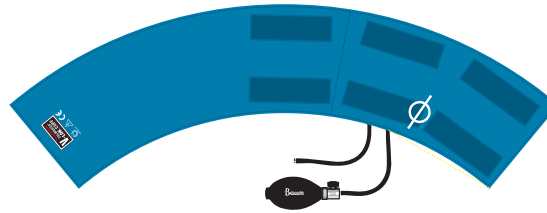

SPECIFICATIONS:

Includes a custom-made, curved Dacron® Polyester Cuff, Natural Rubber Latex or Non-Latex Inflation Bag, Large size Latex or Non-Latex Bulb and Baum Air-Flo® Control Valve.

The BAUM CUFF:

Durable: Baum cuffs are made of tightly woven, heat set, Dacron® polyester fabric that is double stitched at critical points for long wear. Hook and loop fasteners are evenly matched in holding strength to provide the optimal number of open-close cycles. Each cuff is treated with an effective anti-microbial agent.

Comfortable: Baum cuffs and purpose built inflation bags conform to the shape of the limb and have soft edges. Whereas, the use of one piece cuffs may lead to bladder segmentation resulting in uneven compression and discomfort (pinching) of the patient's limb.

THE BAUM BULB & INFLATION BAG:

Latex-Free: Provide a high quality alternative for patients and practitioners who are hypersensitive to natural rubber latex. The inflation bulb and bag are a special blend of Neoprene that is resistant to attack from ozone, acids, bases, alcohols and various hydrocarbons.

W.A. Baum Co. Inc.

SINCE 1916

ORIGINATOR AND MAKER OF ACCURATE
BLOOD PRESSURE APPARATUS EXCLUSIVELY