

# So many choices.

## One standard.

Our meticulous approach to the manufacture of blood pressure instruments has made W. A. Baum Company a leader in the industry. Whether you require wall-mounted, mobile or portable equipment, any of these six instruments will be the right choice for dependable, accurate blood pressure readings.

For more information, visit us online at [www.wabaum.com](http://www.wabaum.com) or call us toll-free at 1-888-281-6061.

Standby®  
Model





Roll-By®  
Mobile  
Aneroid



Desk  
Model



Pocket  
Aneroid

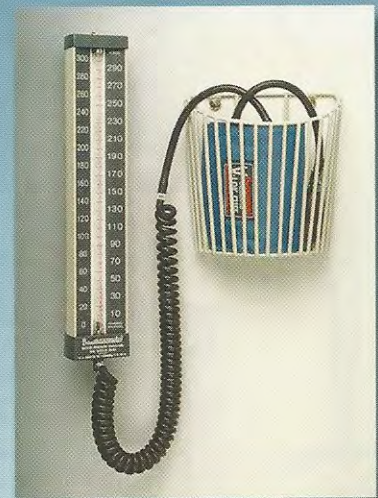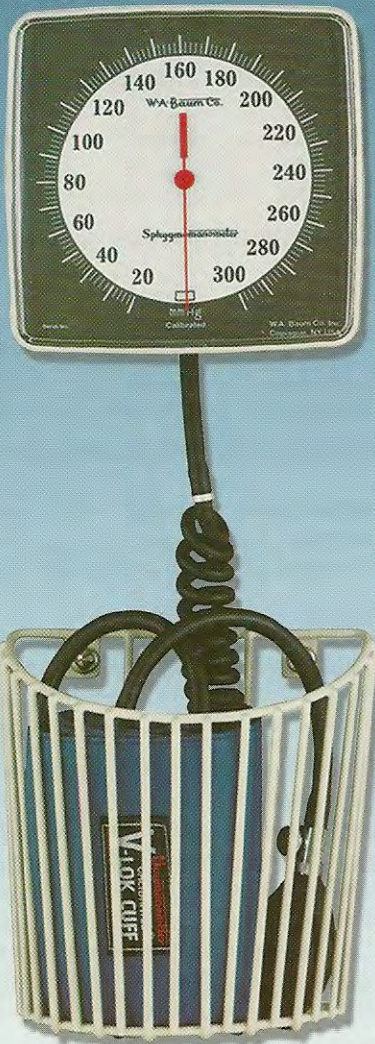

**Baumanometer® Wall Unit 33:**  
Functions on the immutable law of gravity and is guaranteed to be and to remain scientifically accurate.

**Wall Aneroid:** Like all Baum aneroids, it is built to agree with our master instruments, which are traceable to the United States National Institute of Standards and Technology. Each instrument is equipped with our Calibrated® V-Lok® Inflation System.

### W.A. Baum Co. Inc.

SINCE 1916

620 Oak Street, Copiague, NY 11726 USA

[www.wabaum.com](http://www.wabaum.com) • [info@wabaum.com](mailto:info@wabaum.com) • Toll-Free: 1-888-281-6061 • Tel: 631-226-3940 • Fax: 631-226-3969



# Features and Specifications



## The Calibrated® V-Lok® Inflation System:

The Calibrated® V-Lok® system immediately indicates whether the correct size cuff is being used. It consists of a latex bulb, Air-Flo® control valve, latex bag and Calibrated® V-Lok® Cuff. Hook and loop fasteners are used to provide the optimal number of open-close cycles. The center of the bladder is clearly marked on the inner and outer surfaces of the cuff for correct application on either right or left limb. The cuff material is a tightly woven, urethane coated, heat set, Dacron™ polyester fabric. For additional information refer to our catalog, the Calibrated® V-Lok® Cuff data sheet or visit our website at [www.wabaum.com](http://www.wabaum.com).

## Mylar® Clad Calibrated Glass Cartridge Tube:

Since February of 1995, all of our glass cartridge tubes have been enhanced by a new process which involves the application (using heat and pressure) of several layers of very strong and crystal clear Mylar® film which strengthen the tube and maintain its structural integrity even if the inner glass is broken. All of our tubes have a 5mm bore which is 25% greater than required by Federal Specifications; this prevents separation of mercury which, in small bore tubes, cause erroneous readings.

## Wall Unit Ordering Information:

### Cat. No.

| Latex | Non-Latex | Model        | Inflation System         | Calibration |
|-------|-----------|--------------|--------------------------|-------------|
| 0850  | 0850NL    | Wall 33      | Adult Calibrated® V-Lok® | 300 mm Hg   |
| 0950  | 0950NL    | Wall Aneroid | Adult Calibrated® V-Lok® | 300 mm Hg   |

## Mobile Instrument Ordering Information:

### Cat. No.

| Latex | Non-Latex | Model    | Inflation System         | Calibration |
|-------|-----------|----------|--------------------------|-------------|
| 0250  | 0250NL    | Standby® | Adult Calibrated® V-Lok® | 300 mm Hg   |
| 1150  | 1150NL    | Roll-By® | Adult Calibrated® V-Lok® | 300 mm Hg   |
| 1120  | 1120NL    | Roll-By® | Adult Calibrated® V-Lok® | 300 mm Hg   |
| 1160  | 1160NL    | Roll-By® | Adult Calibrated® V-Lok® | 300 mm Hg   |

## Portable Instrument Ordering Information:

### Cat. No.

| Latex | Non-Latex | Model          | Inflation System         | Calibration |
|-------|-----------|----------------|--------------------------|-------------|
| 0320  | 0320NL    | Desk           | Adult Calibrated® V-Lok® | 300 mm Hg   |
| 1050  | 1050NL    | Pocket Aneroid | Adult Calibrated® V-Lok® | 300 mm Hg   |

## Guarantees:

Every Baumanometer® instrument is a precisely calibrated true gravity manometer. It is guaranteed to be and remain scientifically accurate. The entire apparatus is guaranteed against defects in workmanship and material for a period of one year from date of purchase. We guarantee the free replacement of the Calibrated Cartridge Tube if it is ever broken.

Every Baum Aneroid instrument is guaranteed against defects in workmanship and material for a period of one year from date of purchase.

We recommend that every aneroid instrument be checked for accuracy periodically by following our "Testing Instrument Accuracy" procedure\* to determine if it needs adjustment. Adjustments will be made to Baum aneroid instruments without charge when returned, postage prepaid, unless the instrument has been damaged in shipment, tampered with or abused. Our pledge is that our customers will be completely satisfied with our products.

\*Printed copies of "Testing Instrument Accuracy" are available on request or by visiting the FAQ section of our website at: [www.wabaum.com](http://www.wabaum.com).

These guarantees are invalid when other than genuine Baumanometer® parts are used.

Baumanometer, Standby, V-Lok, Extendex, Calibrated, Hi/Lo, Air-Flo and Roll-By are registered trademarks of W.A. Baum Co., Inc.

W.A. Baum Co., Inc.

Tel: 631-226-3940

620 Oak Street

Fax: 631-226-3969

Web: [www.wabaum.com](http://www.wabaum.com)

Copogue, NY 11726 USA

Toll Free: 1-888-281-6061

Email: [info@wabaum.com](mailto:info@wabaum.com)

Since 1916 Originator and Maker of Accurate Blood Pressure Apparatus Exclusively

# W.A. Baum Co. Inc.